


Instructions for use

High-Gold dental Alloy for Porcelain bonding, extra hard, yellow, Palladium free, type 4

Indication

Porcelain-to-metal and conventional application, C&B technique, crown, bridges, milling works.

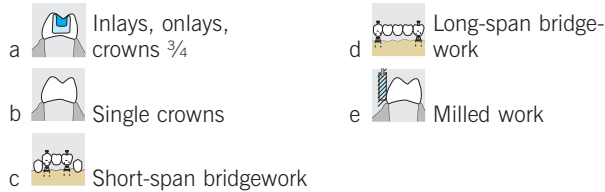

Physical properties

Composition in weight %

Au + Pt group metals	97.85
Au	85.90
Pt	11.70
Ir	0.05
Rh	0.20
In	0.10
Fe	0.05
Mn	0.10
Nb	0.40
Zn	1.50
Density g/cm ³	18.8
Melting range °C	1025–1190
CTE (25–500°C) 10 ⁻⁶ K ⁻¹	14.3
(25–600°C) 10 ⁻⁶ K ⁻¹	14.5
Young's Modulus GPa*	80

Mechanical properties

	1	2	3
Hardness HV5*	185	205	210
0.2% Proof stress, Rp 0.2% MPa*		540	580
Elongation A5%*		6	6

State

1	as cast
2	after firing
3	hardened

Solders before bonding

S.G 1030

Melting range

990–1030°C

Solders after bonding

S.G 750

Melting range

695–750°C

* The values indicated result from measurements obtained under exactly defined conditions. Individual deviations of ± 10% are possible and to be considered as normal.

Preventive measure (contraindication) and important notes

These instructions for use must be precisely followed. It has been suggested in specialized literature that some of its components can, in extremely rare cases, have allergy effects. The choice of the material is the decision of the practitioner based on his knowledge of the sensitivity of the patient.

Corrosion resistance and biocompatibility

V-Gnathos® PLUS has been submitted to the following tests:

Corrosion resistance according to ISO 22674

Cytotoxicity test according to ISO 10993-5

Sensitization test according to ISO 10993-10

Mutagenicity test (AMES) according to 10993-3

The alloy is considered to be highly corrosion resistant and showed neither a cytotoxic nor a mutagenic potential nor did it cause any allergic sensitization.

Launch year 1995

V-Gnathos® Plus corresponds to the standards EN ISO 22674 / EN ISO 9693

V-Gnathos® Plus has been manufactured according to the quality standards ISO 9001 / ISO 13485.

Rx only

The products carry the CE sign.
See packaging for details.

Recommendation

The solder joints surface area should be **10 mm²** to provide sufficient stability to the frame.

Waxing up

The ultimate thickness of the cast metal coping should not be less than **0.3 mm**.

Therefore: Consideration must be given to this fact at the wax-up stage, long span frameworks require structural reinforcement to ensure stability and anticipated solder joints should be of adequate surface area to provide sufficient stability to the frame.

Spruing

Wax sprues of no less than **Ø 3.5 mm** are required. Direct (**Ø 3.5 mm**) and cross bar (**Ø 5 mm**) spruing produce excellent results. Feeder sprues to heavy pontics should be of at least **Ø 4 mm**. Air vents (**Ø 1 mm**) may be used to advantage. Wax patterns should be set outside the thermal centre, i.e. near the casting ring wall and about **5 mm** from the end. For individual copings and small bridges (up to three units), use of the circular sprue provides ideal positioning of the wax patterns and ensure controlled solidification of the frameworks.

Investing

When using steel casting rings always use refractory liner in order to allow free expansion of the investment. All regular or phosphate-bonded (e.g. **Univest®Plus** or **Univest®Rapid**) investments for precious metal alloys may be used. Follow the procedures recommended by the manufacturer.

Preheating

Observe manufacturer's recommendations with regard to setting times, temperature levels etc. On reaching the end temperature a soaking period of **20 to 45 min.** is advisable depending on the size of the cylinder.

Preheating temperature: 750°C

Crucibles

Glaze the crucible before first use with a recommended flux (e.g. Borax/Boric acid). The following crucibles can be used:

Ceramic crucible: Casting temperature 1390°C

Re-melting

When melting down thoroughly cleaned casting buttons, add at least **1/3 of new alloy**. The used copings have to be clean, free of investment and flux residue.

Melting

It is important, when using a torch for melting (for inst. **Meteor type «O»**) that the recommended propane (approx. **0.5 bar** or **7.25 psi**) / oxygen (approx. **1.0 bar** or **14.5 psi**) mixture and pressure are observed. Before melting add a pinch of flux to the alloy. Once the alloy has completely melted, continue heating for a further approx. **5 sec.** before releasing the casting machine arm. When melting by induction or resistance heating, the power must be maintained for an additional approx. **5 sec.**
Flux: boric acid

Cooling

Bench cool the casting ring.

Devesting/Cleaning

Clean by blasting with glass beads or pure aluminium oxide (Al₂O₃) particle size (**50 µm** to **125 µm**) at **1.5 to 2.0 bar** (**21.7 to 29 psi**) pressure.

Pickling

Following casting or soldering, the frame at room temperature may be pickled in hot **Desoxid** for at least **2 min.** – Desoxid I (75%).

Finishing

Trim the framework first preferably using carbide burs and then fine grinding points at low speed. Blasting is done using pure aluminium oxide (approx. **110 µm**).
Recommended grinding points: Cerasiv Brown (5'000–10'000 rpm)
Pressure: 1.0–2.0 bar

Cleaning

Frameworks must be cleaned preceding oxidation either in boiling water for **10 min.** (distilled water is recommended) – in an ultrasonic unit for **5 min.** – or with a steam cleaner.

Gilding (Aurofilm® 2000)

The application of a gilding layer of Aurofilm is possible according to the instructions for use. Aurofilm should be fired in a ceramic oven.

Firing temperature: 920°C–940°C

Oxidation

10 min. under vacuum to 950°C and bench cool. The framework should have a uniform grey-brown appearance. After oxidation frames should be cleaned in hot, diluted (10%) sulphuric acid (H₂SO₄) or in a hot pickling solution of **Desoxid I** for **10 min.** Use this solution just for this alloy. Then thoroughly clean frameworks.

Porcelain application

Porcelain or resin applications are subject to the manufacturers instructions. The alloy is compatible with most known porcelain brands. Consider the CTE of the alloy in conjunction with the cooling recommendation of the porcelain manufacturer. Paste opaques should be pre-dried for approx. **10 min.** before firing at manufacturers recommendation.

Pre-drying temperature: 300°C–400°C.

The porcelain should be fired as recommended by the manufacturer but the first bake should be cooled slowly:
Cooling to 800°C

Excellent results have been obtained with the following porcelains (in alphabetical order), with a normal cooling cycle:
Carat, Ceramco II, Creation &, dSign, Duceram, Duceram Plus, Finesse, MagicCeram, Omega 900, Vintage, Vita M 13

Excellent results have been obtained with the following porcelains (in alphabetical order), with a slow cooling cycle:
Symbio Ceram

Every crown of a bridge must be supported by a supporting stick. An individual firing tray is recommended for long span bridges.
Max. firing temperature: 950°C

Soldering/Laser

We recommend using a propane/oxygen torch (**Meteor Type «L»**) for soldering and a flux like **Fluxor**. During soldering wear dark goggles for protection. The design of the soldering block is a compromise between minimising its thermal mass whilst retaining sufficient strength to avoid its fracturing during soldering. Leave a parallel gap of **0.1–0.2 mm** between surfaces to be soldered and sufficient area to ensure adequate strength of the joint.

Laser welding wire: LW No. 1

Heat treatments

Self-hardened by slow (bench-) cooling of the casting ring
Hardened 450°C/15 min. – let bench-cool

Polishing

The oxide layer may be removed in **Desoxid**. Thorough rubber-wheeling ensures easy polishing. For best prepolishing results, soft brushes, felts and cotton wheels are employed using **Legabril Diamond** diamond-paste. «Whiting chalk» (chalk-powder) mixed with water on soft brushes or cotton wheels may be used for final high polish.